JC-433G – BODY MONITOR

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JC-433G GB

JC-433G – BODY MONITOR

VARIATION OF BODY FAT THROUGHOUT THE DAY

Levels of body hydration, lifestyle and daily habits can influence the detection of the percentage of body fat. A detection performed during the early hours of the morning will calculate a higher percentage of fat because after a long nights sleep, the body tends to be dehydrated. For a more accurate reading and in order to compare different detections it is recommendable to use the product at the same hour of the day and under the same conditions.

Before measuring it is important to remember that:

- The percentage could be lower after having taken a bath
- The percentage could be higher after eating or drinking
- Detections could be irregular during a woman's menstrual cycle
- The percentage of water in our body can be reduced after having been ill or after heavy physical activity, as these conditions make us sweat more than usual.

WHAT IS THE PERCENTAGE OF WATER IN THE BODY?

The percentage of water present in the human body corresponds to the quantity of liquids expressed in a percentage, according to the total weight of the individual, and it is calculated using a particular scientific formula, also used by expert dieticians. It is important to know that the water inside our body performs a fundamental function and is in all tissue and body organs: approximately 50-65% of our body is made of water. Thanks to these scales it is possible to keep the percentage of water under control, which is fundamental for correct functioning of our body. The percentage value of water present in a healthy individual is as follows:

WOMEN from 45-60%

MEN from 50-65%

VARIATION OF THE PERCENTAGE OF WATER THROUGHOUT THE DAY

To obtain a reliable detection it is recommended to measure the percentage of water in the evening (before dinner) as hydration levels are more normal during this hour of the day. In any case, to relate various measurements it is always recommendable to respect established hours. The percentage of water tends to be lower if detected in the morning and this is because the body dehydrates during sleep and above all because during the night liquids tend to concentrate in the torso (whereas during the day thanks to the force of gravity they spread out in a homogeneous manner). Factors that tend to modify the percentage of water are: sleeping, eating and/or drinking, physical activity, having a bath, consumption of medicines, physical conditions and urination.

WHAT IS MUSCLE MASS?

Muscle mass is defined as the complete volume of muscle of the human body. The muscle system is all the tissues that permit the locomotion of the organism itself and the flow of internal organic substances such as blood and food. Control of muscle mass is especially important for people who have problems with weight and fat. Seeing as the body needs muscle to demolish fat cushions, a reduction of weight can give positive lasting results but only if accompanied by the strengthening of muscles.

IMPORTANT WARNINGS AND SAFETY RULES

The values obtained by this body fat monitor are indicative data that **for no reason whatsoever should replace a doctor's opinion**. Do not forget: self-measuring means control, not diagnosis or treatment. Unusual values should always be discussed with your doctor.

Under no circumstances should you modify the dosage of any medicine prescribed by your doctor.

The scales should not be used by:

- People with pacemakers or any other internal device as the weak electric signal could in fact interfere with correct functioning of the device.
- People with oedema, osteoporosis or fever
- People in dialysis
- If taking medicine for cardiovascular problems
- Pregnant women
- Athletes that practise more than 10 hours of intensive training per week
- People who have a resting pulse frequency that is equal to or below 60 beats per minute.

Only use the scales for the intentions indicated in this instruction manual and exclusively for domestic use. Any other type of use is to be considered improper and therefore dangerous. Check that the appliance is whole and is not damaged before use. If in doubt, do not use the product and consult a specialised after sales centre. To not attempt to repair the product yourself, always go to a specialised after sales centre. To avoid any interference it is recommendable to use the product away from appliances with a strong magnetic field.

Never place the scales on wet or slippery surfaces. Never immerse the product in water and never put it away in a damp place. Always keep the product away from sources of heat, direct sunlight, cutting objects or similar objects.

When the appliance is not used for a long time, it is recommendable to remove the battery as leaking of fluids could damage the product. Also, seeing that the quality of batteries varies between models, the warranty does not cover damage caused by any possible leaking of fluids. Never leave the product unattended when children are present and always put it away out of their reach.

FEATURES

Capacity: 180kg

Division: 0.1kg / 0.2lb

Weight unit: kg / lb

Minimal weight: 10kg

Function key: ◀, SET, ▶, UNIT

Display: Blue LCD with white backlight

Battery: 4xAAA

User No.: 4

Age range: 6~100

Stature range: 80~220cm / 2ft 7.5inch~7ft 2.5inch

Athlete measuring mode

User identification: Identify user automatically according to the last measuring weight


Body fat, water, muscle and bone measuring function

Fat grade indication

Auto on or press UNIT or touch SET to switch on

Auto off function

LCD display

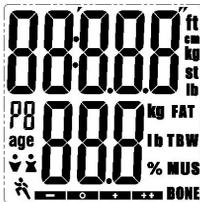Symbol	Description
	male
	female
	athlete
	underfat
	healthy
	overfat
	obese

Function indication

- ◀: Touch "◀" to decrease the value in setting mode.
- SET: Touch "SET" to confirm the value in setting mode.
- ▶: Touch "▶" to increase the value in setting mode.
- UNIT: Press "UNIT" to select desired unit in weighing mode.

Operation instruction


I. Select unit before use

a) Open the battery door at the back of the scale.

b) Put the batteries in and ensure +/- terminals are correct. The LCD will display .

c) Close the battery door and wait until is shown on the LCD.

d) Press "UNIT" to select the desired unit.

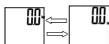Note: If there is not any operation in weighing mode, the scale will turn off automatically after about 10s.

II. Input personal information before use (the unit kg selected)

a) Touch "SET" to switch on or touch "SET" in weighing mode. The user No. Pn (n can be 1~8) will flash.

b) Touch "◀" or "▶" to select the desired user No. and then touch "SET" to confirm to enter setting mode.

c) Touch "◀" or "▶" to select gender or athlete mode and then touch "SET" to confirm.

d) Touch "◀" or "▶" to set up the value of stature and then touch "SET" to confirm.

e) Touch "◀" or "▶" to set up the value of age and touch "SET" to confirm. The data would be saved successfully, and then you may start measuring after

is shown on the LCD according to the measuring procedure below.

Note: If there is not any operation in setting mode, the scale will turn back to weighing mode after about 20s.

III. Measuring(select user manually)

a) Please put the scale on flat hard ground.

b) Touch "SET" to switch on and use "◀" or "▶" to select your No..

c) Please wait until the LCD display "0.0".


d) Stand on the platform with bare feet.

e) Stand still and the weight data will be locked automatically when stable.

f) The measuring will start automatically with the "o" signal moving from left to right.

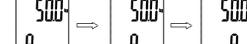g) After finishing measuring, the results will be shown in cycles for three times and then the scale turn off automatically.

IV. Measuring(user identification automatically)


a) Please stand on the platform directly with bare feet to switch on.

b) Stand still and the weight data will be locked automatically when stable

c) The measuring will start automatically with the "o" signal moving from left to right.

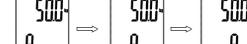d) If there is only one user whose last weight data is close to the current one, it will be selected automatically. The results will be shown in cycles for three times and then the scale turn off.

e) If there are more than one users whose last weight datas are close to the current one, they would be shown on the LCD in cycles after calculation. For example

f) Please use "◀" or "▶" to select the user. Touch "◀" to select "P1" and "▶" for "P2".

Note: If there is not any selecting operation, the scale will turn off after about 40s.

g) After finishing selection, the results will be shown in cycles for three times and then the scale will turn off.

Note: Before using automatic user identification measuring method, please select your No. to measure once at least, or failed identification.

Indication message

a) Low battery indication. When "Lo" is shown, the scale will turn off after a few seconds and please change the batteries.

b) Overload indication. Please remove the weight on the scale to protect the loadcells.

c) Measuring error indication. The weight data is shown only and the scale will turn off after a few seconds. Please stand still with bare feet while measuring again.

d) Failed identification indication. The weight data is shown only and the scale will turn off after a few seconds. Please select your No. to measure once again.

Reminder

- a) Please do not put the scale on carpet.
- b) Please stand on the platform with bare feet to making the correct measuring.
- c) If the scale is not to be used for a long period, please remove the batteries.
- d) Clean the scale with a slightly damp cloth, DO NOT immerse the scale in water or use chemical/abrasive cleaning agents.

Body fat, water, muscle range

a) The body fat ratio (%): 3%-50%

Age	Female				Male			
	Underfat	Healthy	Overfat	Obese	Underfat	Healthy	Overfat	Obese
<39	<21	21.1-33	33.1-39	>39	<8	8.1-20	20.1-25	>25
40-59	<23	23.1-34	34.1-40	>40	<11	11.1-22	22.1-28	>28
>59	<24	24.1-36	36.1-42	>42	<13	13.1-25	25.1-30	>30

b) The body water ratio(%): 25%-75%

Age	Female			Male		
	Low	Normal	High	Low	Normal	High
10-15	<57	57.1-67	>67	<58	58.1-72	>72
16-30	<47	47.1-57	>57	<53	53.1-67	>67
31-60	<42	42.1-52	>52	<47	47.1-61	>61
61-80	<37	37.1-47	>47	<42	42.1-56	>56

c) The body muscle ratio(%): 25%-75%

Age	Female	Male
	Healthy	Healthy
6-100	>34	>40

CLEANING AND MAINTENANCE

Use a soft cloth to clean the platform. It is also possible to use a slightly dampened cloth. Do not use chemical or abrasive substances. After cleaning the platform make sure it is completely dry before carrying out a measurement. If the surface is wet or damp it could become slippery. Always position the scales on a stable and flat surface. Avoid rugs and similar objects.

Avoid jumping on the scales, dropping them or subjecting them to particular knocks. As well as for precision, the internal mechanisms are in fact very delicate, and could easily be damaged. Avoid exposing the product to electromagnetic waves that could cause the loss of data in the memory.

The seller is not responsible for loss of data memorised in the scales.

Note: All indications and diagrams are based on the latest information available when the manual was printed and may be subject to change.

WARNINGS

These scales should not to be used to weigh objects or substances in commercial transactions, to produce medicines, to calculate tolls, tariffs, taxes, premiums, fines, remunerations, indemnities or fees of a similar kind determined by weight. These BFW scales are to be used for keeping one's own personal body weight,

and percentages of fat and water, under control and are not to be used for making diagnoses or for medical treatment. Values that may seem unusual or outside the norm should always be discussed with your doctor.

In order to avoid danger to the user of this device, the doses of any drug or therapy prescribed by one's own doctor should in no circumstances be modified on the basis of data calculated by this device.

Specifications and designs are based on the latest information available at the time of printing and subject to change without notice.

The product is made conforming with all the applicable European regulations.

DISPOSAL

The device (including its removable parts and accessories) must not be disposed of together with municipal waste at the end of its life, but in compliance with European Directive 2002/96/EC. Since it must be handled separately from household waste, either carry it to a separately-collected waste disposal centre for electrical and electronic appliances or give it back to the retailer on purchasing a new device with the same purpose. Any infringement will be severely prosecuted.

The batteries used in this device must be disposed of in the special bins at the end of their life.

NOTE: Specifications and designs are based on the latest information available at the time of printing and subject to change without notice.

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GB: WARRANTY CARD

Joycare S.p.A., a company with a sole shareholder, offers a period of 3 year warranty after the date of purchase. The warranty does not cover damage caused by accidental breakage, transportation, incorrect maintenance or cleaning, misuse, tampering or repairs carried out by non authorised third parties. The following is also excluded: damage caused by incorrect product installation, and damage caused by wear and tear, obviously including parts subject to wear and tear. In order to obtain the replacement/repair of products under warranty in cases other than those mentioned above, but which are instead caused by manufacturing defects, please contact your local importer/distributor. Goods that are not held to be in good condition because of the aforesaid causes must be accompanied by this certificate, duly filled-in, as well as a receipt for the goods (invoices and transport documents) and purchasing documentation.

Please keep the original packaging of the product.

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JOYCARE S.p.A. socio unico